



IMPLANT ONE IMPRESSION POST, OPEN TRAY *(Impression Pin, Impression Coping)*

Used to make open tray impressions.


Made from 6AL4V E.L.I. Titanium.

Part Numbers


300 series	400 series	500 series
IPO-0300-00	IPO-0400-00	IPO-0500-00



USAGE

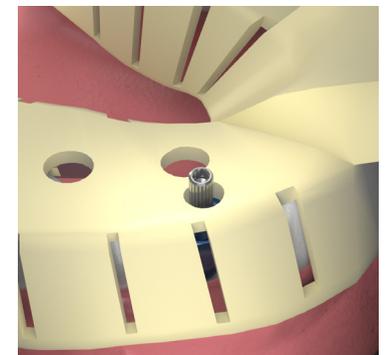

Customize the impression tray so that the center pin sticks out through it.

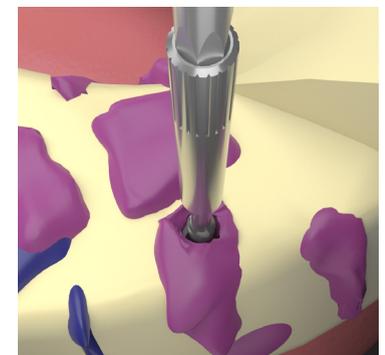

After impression has cured, remove the center pin with the .050 hex tool. The body of the impression post will remain in the impression material when the tray is removed.

Center pin secures post to implant, and is removed after the impression has cured.



Post body remains lodged in the impression material.



Flat face is parallel with hex flat.



Specifications

Material: 6AL4V E.L.I. Titanium
Profile Diameter: 3.2 mm – 4.7 mm

